



MEDILAP


Quality assures reliability & guarantee...



MEDILAP 400 MAESTRO

HIGH PERFORMANCE ELECTRO SURGICAL UNIT

MEDILAP Electro Surgical Unit Features:-

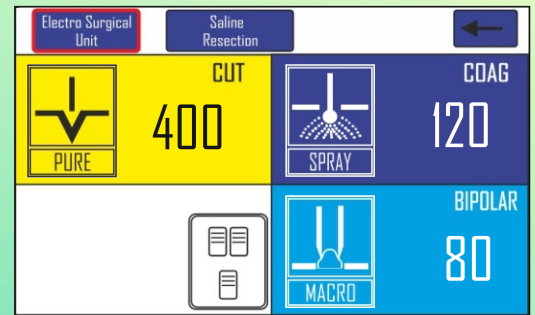

Monopolar Modes with two Individual output ports-

Two independent power activation of monopolar output, hand switch operation, and foot switch operation. The auto-monitoring facility displays different error codes if any error occurs during use. Randomized form of spray mode.



Monopolar Modes

Up to 400 W Monopolar Power for cutting and 350 W for coagulation with 4 independent Cut/Coag modes, with precise power-adjustable using a digital interface.

Monopolar Cut (Pure Cut)

for a smooth cut **with minimal** eschar formation. This so-called “pure cut” comes closest to the classic scalpel cut, thus providing the best basis for effective wound healing.

Monopolar Cut 2 LowCut)

for a pure cut with low operating frequency and power.

Monopolar Cut 3 (Blend Cut)

Blend II or Super Blend is a cutting current with even **eschar formation compared to** Blend Cut. As the wound edges are simultaneously coagulated during the cutting process, bleeding from dissected capillary vessels is stopped at once and extended hemorrhages are prevented.

ENDO Mode


This is a time-controlled cutting mode offering a fractionated and therefore controlled cut intended for special applications, e.g. for endoscopic papillotomy and endoscopic polypectomy.

Monopolar Soft (Deep)Coagulation


Coagulation current with depth effect (electrode in direct contact with the tissue).

Monopolar Forced(Speedy) Coagulation

Forced coagulation current enabling fast coagulation. An ideal current for dissecting tasks because it provides effective cutting capability despite its dominant coagulation feature.

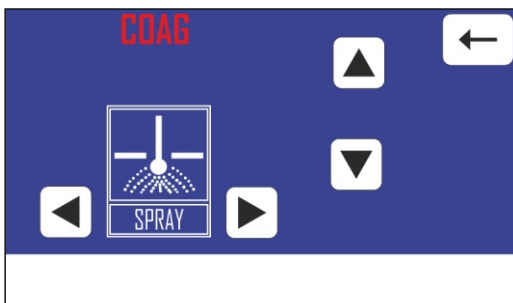


MEDILAP Electro Surgical Unit
Features:-

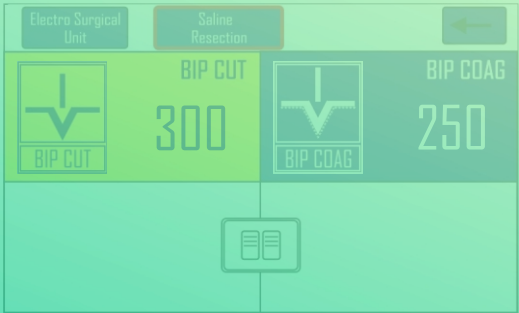


Bipolar Modes - Consist of Auto/Manual mode

4 different Bipolar modes give users additional flexibility and precision while performing any procedure. Separate bipolar generator with foot switch operation or hand (auto) operation, eliminates the need for foot switch while using auto mode.

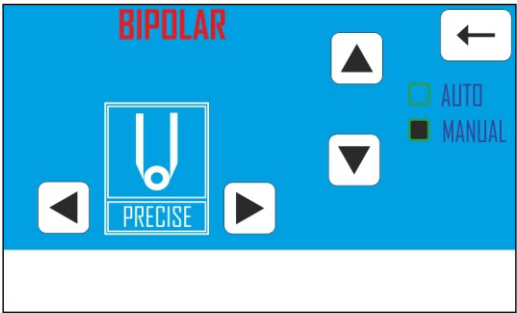

Bipolar Mode with Separate Output Port

Up to 120 W Bipolar Output, consist of 4 different modes for bipolar cut and coagulation, providing the user easiness to adjust the power as per the requirement.

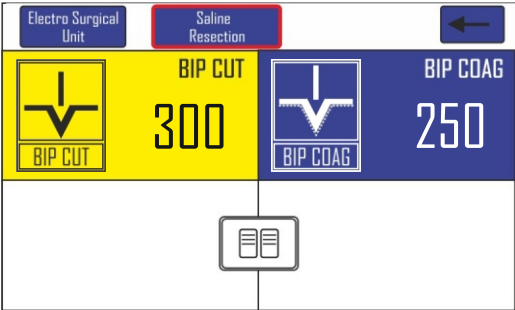

BIPOLAR(Cut)

The Bipolar modes used to cut tissues with the help of bipolar current. Special bipolar forceps/tools are used for this purpose. An increased level of haemostasis is more effective in terms of controlling bleeding, but it also affects the tissue thermally more.

BIPOLAR(Coag)

In Bipolar coag mode, the current flows between the electrode tips and there is no need to apply the neutral electrode. A typical application is closing medium size blood vessels.



Patient plate fault monitoring system.

Excellent multi functional high frequency electro surgical device

User convenient facilities and settings with large visual overview of touchscreen display which makes this device versatile.

Auto monitoring facility, indicates different error codes if any error comes during use.

Randomized form of spray mode.

Advance Cutting Mode Pure Cut with tissue sensing technology

Two independent power activation of mono polar output, hand switch operation and foot switch operation.

A tissue sensing technology used for all applications to deliver user consistent performance(real time feedback system).

System for neutral plate safety by continuous monitoring of contact quality and connection.

Separate bipolar generator with foot switch operation or hand (auto) operation, eliminates the need of foot switch while using auto mode.

Shared Monopolar Cut & Coag output shared between two activated monopolar instruments at the same time.

Generator is compatible with laparoscopic monopolar and bipolar instruments.

Auto Bipolar Mode with timer up to 7sec

INPUT POWER SUPPLY

Input Voltage 180-270V Ac
Current Fuse 6 Amps

SAFETY FEATURES

Unit Construction According IEC-60601-2 & 60601-1
Low Frequency Leakage Currents According to IEC 60601-1
High Frequency Leakage Currents According to IEC 60601-2-2

Protection I
Unit Type CF
Output All Floating Output
Power Efficiency > 98%
Nominal Frequency 350 kHz
High Frequency <150μA RMS in all output
Low Frequency <10μA max.

MONOPOLAR CUTTING

Pure Cut 400W for 500 Ω Resistance
Low Cut 300W for 500 Ω Resistance
Blend Cut 350W for 500 Ω Resistance
Endo Cut 250W for 500 Ω Resistance

MONOPOLAR COAGULATION

Soft 200W for 500 Ω Resistance
Dessicate 150W for 500 Ω Resistance
Force 175W for 500 Ω Resistance
Spray 120W for 500 Ω Resistance

BIPOLAR

STD 150W for 100 Ω Resistance
Macro 125W for 100 Ω Resistance
Micro 100W for 300 Ω Resistance
Bi-cut 150W for 100 Ω Resistance

BIPOLAR SALINE RESECTION

BI CUT 300W for 100 Ω Resistance
BI COAG 250W for 100 Ω Resistance

OTHER FEATURES

Constant control Real time Monitoring
Volume control 4 Levels
Multi Switch Separate Monopolar & Bipolar Foot switch
Output Port 2 Monopolar & 1 Bipolar

WEIGHT & DIMENSION

Length 390 mm.
width 340 mm.
Height 130 mm.
Weight 6.5 Kg.

POWER SUPPLY

Input voltage 230v 50Hz +/-10% or 130V 60Hz +/-10%
Related power output 500W